## **Addition to the Assembly instructions!**

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The Swedish Back Care System™

## **Back-A-Traction®**

**ONLY pertains to**  
**Model CN-AH6 Heavy Duty**

### **Important!**

Before you use Back-A-Traction, model CN-AH6 Heavy Duty, the enclosed “Triangle Enforcements” are to be attached on each side of the frame (A), as shown in the picture below.

Hang the “Triangle Enforcements” over the sides of the stands – see below pictures 1 & 2.

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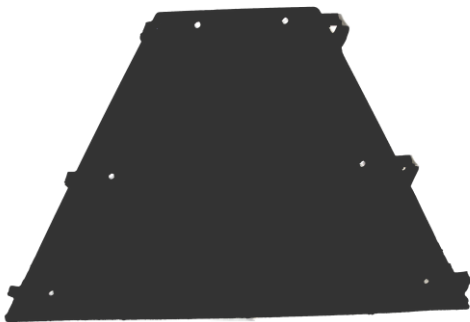
2)

3)

## **Enclosed Lumbar Support Pillows with Model CN-AH6 Heavy Duty**

With Mastercare Back-A-Traction, Model CN-AH6 Heavy Duty, there are 3 Lumbar Support Pillows in different heights included (A), as shown in the picture below.

There are also 3 different positions in the back-rest to attach one of the Lumbar Pillows – to make it fit for best support and comfort (1 and 2).

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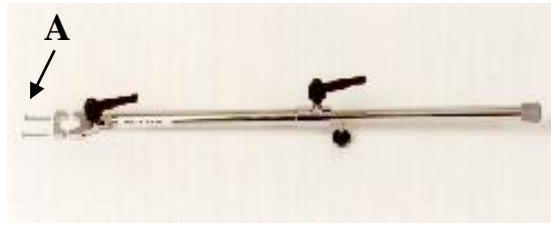
**B**

Place the backrest in horizontal position. Undo the screws (A).  
Place the Support Leg in the middle of the tubular frame, underneath  
the cut-out for the face in the backrest (B), see above pictures.

Fasten the screws in place, the Support Leg will hang straight down.  
The Support Leg can hang loose ready for use (C)...or can be tightened,  
with the handle, in order to hang underneath the backrest (D) when not in use  
(see pictures below).

**Always use the handle (1) – never pull the Support Leg itself !**

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**D**

**1**

# Operating Instructions: Support Leg, SL1

## Accessory to Mastercare Back-A-Traction, Model A

Mainly for Professionals as support during use for massage, manipulation and/or adjustments. Can also be used as support for heavy people (over 220lbs./100 kg, but max 440lbs./200kg), as well as to reduce angle of inclination.

Please note; for Professional treatments in Mastercare Back-A-Traction the Support Leg must be used in order for the warranty to be valid.

### **To properly use the Support Leg you must use the handle.** **It is the handle that fastens or loosens the Support Leg.**

#### **The upper handle - see picture above (1):**

If you pull the handle straight out (with thumb against button) and keep it there, you can rotate it to any position you want...when the handle is back in place you can lock or loosen the Support Leg.

#### **The bottom handle - see picture above (2):**

If you loosen this handle you can lengthen the Support Leg and thus reduce the inclination of Mastercare; for example to only 10 or 5 degree inclination or to horizontal position.

Please note; when you use the reduced inclinations, the locking-devices will not function, but with the help of the bodyweight you will stay in place.

*If either handles are pointing out and are "in the way"...simply pull the handle straight out and keep it there and rotate the handle to desired position.*

**Always use the handle –  
Never pull the Support Leg itself!**

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**Swedish Backcare System, Inc.**